

FIG1 THERMAGDELUXHP HIGH POWER INSTRUCTIONS

The buyer may want to consult with a doctor or qualified health practitioner who is familiar with alternative and complementary medicine before using this product for any particular ailment. We do not claim to be doctors in medicine and are unable to validate or substantiate the claims for cures as stated or implied. We cannot assume any liability for claims, uses, and experimentation. We are presenting this data that we believe is real and authentic as it is from qualified researchers in these fields. Device is not FDA approved.

IF WE HAVE MADE ANY ERRORS BY COMMISSION OR OMISSION, PLEASE LET US KNOW.



See FIG 2, 3 for control identification and functions

Uses Conditioned Magnetic Pulses in Health Research for Increasing Blood Flow and Cell Regeneration, Growth Enhancement Related to Sexual Dysfunctions, Treatment for Many Symptoms of Pain Relief


NOTE Decals are included with the unit and can be put on by the user if desired. Many customers have requested not to label the controls. We therefore leave this up to the purchaser. You may use figure 2 for identification and location.

NOTE Unit is shipped with our basic coil that can be used for just about all function to some degree. Effect will vary due to shape and size. Other coils are available but it is best to tell us what your use is and we can make one up. We have several other standard models as shown. You may want to consult a professional in the use of pulsed magnetic therapy or your personal physician

1. Familiarize yourself with the control descriptions on page 2, 3 of these instructions
2. Read below data on a more detailed explanation of controls and features

FRONT PANEL CONTROL FUNCTIONS

TMPHP1PAGE2

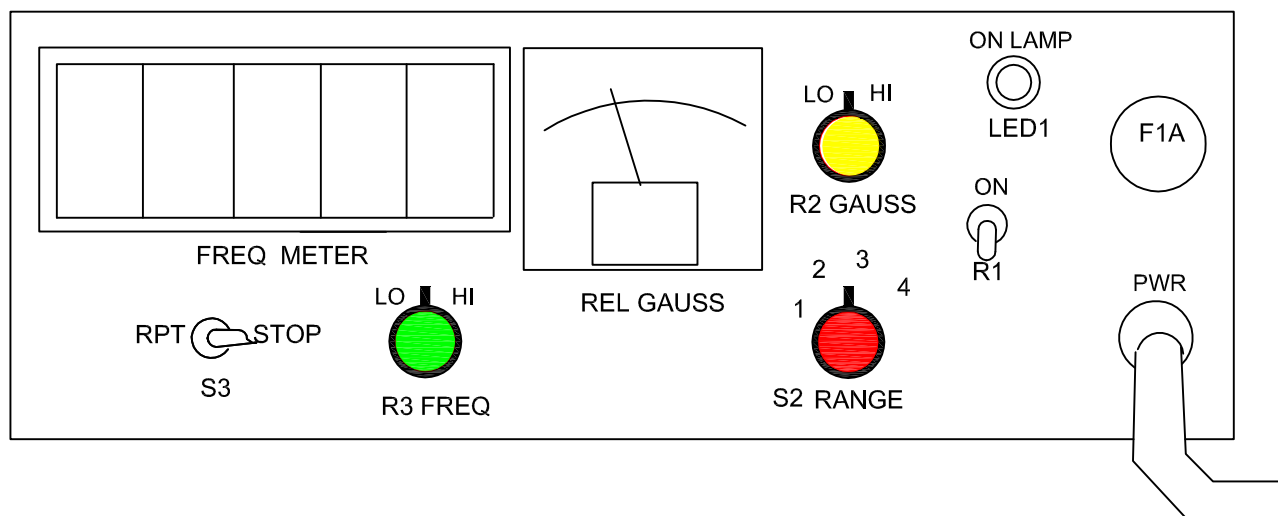

S1 PWR SWITCH turns unit on in up position.

S2 RANGE SWITCH pos 1 selects "10 sec" and pos 2 selects "10 min" for sweep time of sleep and relaxation states. This *sleep cycle* function automatically causes the frequency to start at 100 hz and drop down to 1 hz over the selected period of time of 10 sec or 10 min. **Position 3 selects a freq range of 1-100 hz . Position 4 selects a freq range 50-1500**

S3 CYCLE SWITCH selects this function to keep repeating or stopping after the first sweep of the *sleep cycle*. The freq change is read by the freq meter.

R3 FREQ CONTROL is a 10 turn precision control pot for accurately adjusting frequency ranges as selected by **S2 RANGE SWITCH** pos 3 or 4 and displayed on the **FREQUENCY METER**

R2 GAUSS supplies voltage to energy cap producing $Joules=1/2CVexp2$ L of head=50UH, Capacity=4.4 uFD, PEAK VOLT =450, Peak Current=120 AMPS, Turns on coil=60, Gauss pk=>6000 PLUS!!

LED1RED- Power "on" when S1 power switch is on pos 2

J1-4 banana jacks for hooking in two coils on rear panel

F1- 1 amp fuse holder

PWR 3 conductor #18 power cord with green wire securely attached to chassis. Cord must be plugged in to a 3 wire system for obvious electrical safety


FIG 3 REAR PANEL AND ADDITIONAL NOTES

The meter on the front panel reads the voltage on the energy storage capacitor . This capacitor must charge to the selected voltage for the required energy to produce the magnetic Gauss. This charging function must occur during the period between the pulses. The discharge is thru a matching transformer and then to the terminals where the coil head is plugged in.

There are two frequency ranges on the unit controlled by S2 range switch; POS3 selects 1-100 Hz while POS4 selects 50-2000Hz.

1. In POS 3 you will be able to adjust for max Gauss through the entire frequency range of 1-100hz
2. In POS 4 you must adjust Gauss at $\frac{1}{2}$ meter scale or unit will over heat. Note above 1500 Hz that Gauss meter will start to drop but unit is still functional at those higher frequencies.

You may experience some heating in the coil. If this gets to point of being hot, you turn down R2 gauss control until coil reaches a temperature of comfortably warm

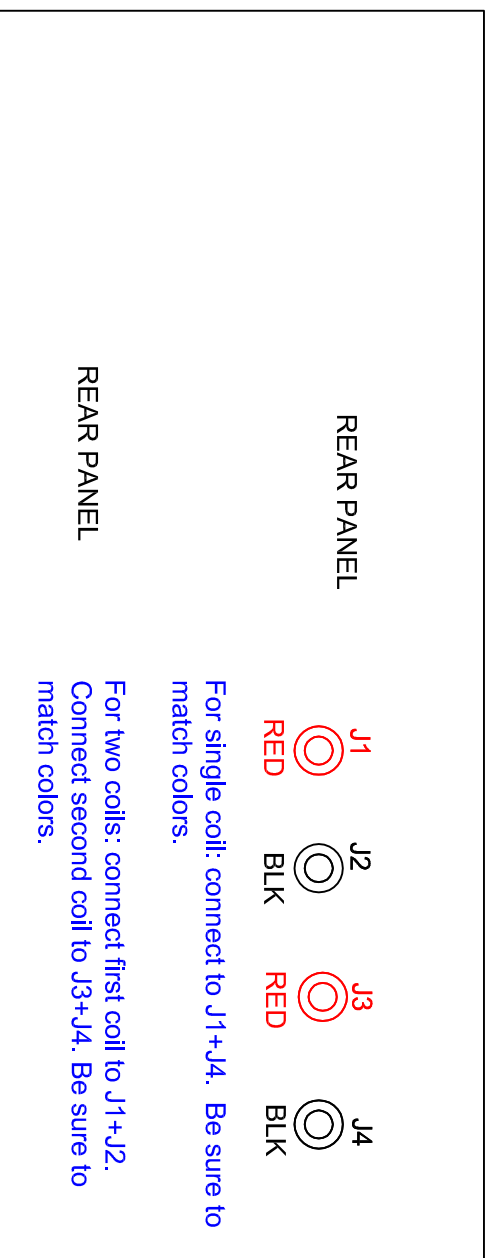

OPTIONAL COILS for use with the THERMAG DELUX HIGH POWER VERSION coils id marked for north pole

Please note that there can many possibilities of coil sizes and configurations dependent on the target treatment area and location.

It is suggested that you send us a sketch indicates size and shape. **CAUTION** use for head or neck applications you must consult a qualified practitioner or doctor.

See below available coils. All have proper mating plugs:

1. Standard 4 inch 6 kilo-gauss flat Archimedean un-potted. Included with unit. extra coil for sandwiching or for second user...39.50.
2. Solednoid coils for wrists, hand, fingers and erectile dysfunctions and size enhancement
- 3 Flexible flat 12" diameter for sleep function inserts under pillow.
4. Room area coils up to 10 x 10' for study of plants, cages or do entire bedroom etc
5. Large solenoid from 2' diameter to 4 ' height.
6. SEPARABLE coil for neck. **caution**
7. Coils for heads. **caution**.

We will be introducing our super high powered pulser of 1-10 pps with energy of 100 joules per pulse. Intended for behavior modification studies

and full size magnetic cage to accommodate a full grown adult. These will only be available to qualified practitioners

The THERMAGDELUX is an all new 115 vac/ portable or battery powered Magnetic Pulser. An optional cigarette plug is available for vehicle use. Those familiar with our lower cost THERMAG10 will appreciate the many added features to the THERMAGDELUX. The THERMAGDELUX MAGNETIC PULSE GENERATOR produces precise selectable low or medium power magnetic pulses. These pulse repetition rates can be selected between 1Hz to 2000 Hz per second via a 10 turn micro-pot. A digital frequency meter displays the selected repetition rate. *A special circuit provides a sweeping action where the pulses automatically change over a selectable period of 10 minutes or 10 seconds. These pulses start at 100Hz and vary down to 1 Hz and are intended for sleep and relaxation enhancement.

****Unique "automatic sleep and relaxation mode" sweeps you through the alpha, beta, theta, delta frequencies at either a 10 or 1 minute full cycle period. You may choose this Auto Sleep Sweep™ feature to auto reset and repeat or just self-terminate after you hit the land of "nod"***

The basic coil that comes with the unit is 4 inches in diameter and contains a Velcro strap for securing to most target points of the body. It comes with a detachable coiled extension cable that can be attached to other optional magnetic coils for other functions

Magnetic Pulser is designed so the North or South Pole can be applied to the body only. The coil strength and polarity is easily measured with our MX3511 MAGNETIC STRENGTH and POLARITY METER found on our web site www.amazing1.com

An interesting example is when an aluminum beer or pop can is held between your ear and the Pulser coil. You now can hear a strong distinct pulsing as the can itself becomes the opposing polarity. If you hold the coil in front of a TV you can see lines forming on the screen. "DOES NOT WORK ON PLASMA OR LCDS" On some TV's the lines are visible when it's held as far away as 3 inches. This magnetic pulsing action is more powerful, and effective, than traditional magnetic therapy.

The pulsing coil takes magnetic therapy to a whole new dimensional level, with exciting break through results and benefits. Permanent magnets, no matter how strong will not produce the same results as the pulsed fields of this device, which produces induced back-emf currents. You must have a high-intensity, time-varying, magnetic impulse, not just a magnetic field. Magnets have their uses, but they work in different ways with different results. Equally important, is the fact that we have "electric-powered" brains: all our thoughts and perceptions consist of complex networks of electrical signals and electromagnetic fields that pulse and sweep throughout the brain. So it then makes sense that harmonic electrical revitalization of the brain can harmonically influence your mental state and positively alter mental effectiveness.

Currently fixed magnets are sold for health applications and are accepted by the medical profession. These magnets only produce a very weak stationary field and cannot reach to the depths in tissues to be as effective

as the pulsed devices. Many users of pulsed magnetic devices using this therapy have claimed a remarkable healing effect on many ailments and afflictions. This data is both from professionals to individual users.

Magnetic pulser's have been registered as a Canadian medical device.

The THERMAGDELUX Pulser as sold is non-invasive yet is not FDA tested. We now offer the device as a fully functional product built that is suitable for experimental use.

Other Benefits

Externally applied magnetic pulses to the lymphatic system, the spleen, kidney & liver, help neutralize germinating, latent and incubating parasites

of all types, helping to block re-infection. This speeds up the elimination of disease, restores the immune system, and supports detoxification.

The movement of the lymphatic system is essential in purifying, detoxifying and regenerating the body, supporting the immune system, and maintaining health. Normally the lymph is pumped by the movement of our body's electromagnetic field with vigorous exercise and physical activity. However, a clogged, sluggish, or weak lymphatic system prevents the body from circulating vital fluids and eliminating toxic wastes, thus weakening the immune system. It makes us vulnerable to infections and diseases. In order to be healthy, it is essential to keep the energy balanced and fluids moving so that the body's natural intelligence may operate in its full healing capacity. In addition, each cell must be enlivened with its own unique

energetic frequencies and harmonic energy-state, and be harmonically connected to the life-force energy throughout the rest of the body.

Pulsed Electromagnetic Fields influence cell behavior by inducing electrical changes around and within the cell. Improved blood supply increases the oxygen pressure, activating and regenerating cells. Improved calcium transport increases absorption of calcium in bones and improves the quality of cartilage in joints, decreasing pain dramatically. Acute and even chronic pain may disappear completely.

This device uses complex energy pulses of magnetic waves to stimulate certain body functions. It accelerates production of vitally important hormones often providing miraculous effects. Magnetic Pulsing aids in human growth hormone production and neurotransmitter production. The results are remarkable, increasing vitality, sexuality, **male penile growth in size**, learning and reduction of memory loss with even reports of increased psychic ability for some people. It revitalizes and rejuvenates. As we get older, hormone production drops off and is why we age.

Users are experiencing faster healing of injuries, including bone fractures, carpal-tunnel and drug free arthritis pain relief. The pulsed magnetic field stimulates blood flow and cellular respiration in the area applied. There are reports of migraine headaches ceasing after magnetic therapy is applied. This pulser can be worn, held or wrapped to any part of the body with your hand or held in place with a

tension bandage. You can use your pulser while driving and on the go, rather than having to stop and apply it when sitting or lying down as is necessary when using a 110 volt only standard Magnetic Pulser.

The following list of pulse frequencies were provided from the below web sites and are not a result of our research. They are claimed to be those used for great healing--and experimental—tools in pulsed magnetic force therapy. These are by no means all of them as research is finding more

Normalize Adrenal Function 1335 Hz

Normalize Pituitary Function 635 Hz

Stimulate Increased / Normalized HGH Production [Pituitary] 1725 Hz, 645 Hz, 1342 Hz

Stimulate Normal Pineal Function 480 Hz

Stimulate Normalized Hypothalamus Function 1534 Hz, 1413 Hz, 1351 Hz

Normalize Endocrine System Function 1537 Hz

Stimulate / Normalize Immune System Function 835 Hz

Stimulate Normal Colon Function 635 Hz

Stimulate Normal Thyroid Function 763 Hz

Normalize Progesterone Levels ; (in sequence) 763 Hz, 1446 Hz, 1443 Hz, 763 Hz

Normalize Estrogen Production Levels : male & female 1351 Hz

Normalize Testosterone Production Levels: Male 1444 Hz

Normalize Testosterone Production Levels: Female 1445 Hz

Stimulate Normal Pancreas Function 654 Hz

Stimulate Normal Liver Function 751 Hz

Stimulate Normal Kidney Function 625 Hz

Stimulate Normal Heart Function 696 Hz

Normalize Blood Pressure 15 Hz

Stimulate Normal Nervous System Function 764 Hz

Stimulate Normal Lymph System Function 676 Hz

Stimulate Increased Lymph System Circulation 15 Hz (15.2 Hz)

Stimulate Normalized Blood Circulation 337 Hz

Stimulate Increased Blood Flow / Circulation 17 Hz

Normalize Red Blood Cell Production 1524 Hz

Normalize White Blood Cell Production 1434 Hz

Normalize Hemoglobin Production 2452 Hz

Stimulate the Reinforcement of DNA Integrity 528 Hz

Stimulate the Reinforcement of RNA Integrity 637 Hz

Stimulate Clarity of Thought / Mental Function 35 Hz

Stimulate the Stabilization Of Emotional States 15 Hz

Stimulate The Clearing of Emotional Trauma / Energy Blocks 15 Hz

Stimulate The Balancing Of Spiritual Well-Being 1565 Hz

Reduce Chemical Sensitivity 443 Hz

Reduce Electrical Sensitivity 657 Hz

Stimulate the Normalization of Calcium Metabolism 326 Hz

Stimulate Repair / Healing Of Nerve Damage 578 Hz , 764 Hz, 657 Hz, (also try 2.0 Hz, 657 Hz,

Stimulate The Healing Of Bones 7.0 Hz

Stimulate the Healing Of Ligaments 9.7 Hz

Stimulate The Healing Of Muscles 13.5 Hz

Stimulate the healing of Capillaries 15.2 Hz

Reduce Swelling of Herniated Disc 25.4 Hz; 326 Hz, 15.2 Hz

Reduce Excess Fluid Retention in Joints & Tissues 15.2 Hz, 24.3 Hz

Reduce General Back Pain (Fibromyalgia) 326 Hz to 328 Hz

Stimulate Indocin Production (stroke recovery, etc.) 2642 Hz

Accelerate the healing & clearing of scarring 5.9 Hz

Accelerate Healing of injuries and Surgeries; Promote Regeneration 47 Hz

Reduce effects of Diabetic Neuropathy 73 Hz

See more at www.elixa.com/light/nogier.htmlFrench Neurologist has a tested list

DREAM ENHANCING (EXPERIMENTAL)

We have received feedback from some of our users” that dreams appear to be programmable by certain combinations of the settings when used in the “*auto sweep mode” **This feature is explained below.**

This could be a new concept in magnetic pulse technology. Programmed pulses of time variant contra-polar resonant magnetic energy are exposed to the brain area resulting in increased neuron activity when in a deep sleep mode. The results appear to produce enhanced dreams of intense reality.

Our research is currently attempting to gather more of this data and analyze it. We then will try and see if there is a correlation between the dreams and the settings when used in this “auto sweep mode”.

It would be great if one could select the dreamed material and actually program an erotic, flying or subject of choice. This research is showing signs of dream repeatability but we are not claiming this feature at this time. Once we succeed, we will include in our product advertising.

Special precautions are the possibility of intense dream resulting in very realistic nightmares. Some of our test subjects have experienced this phenomena, some more than others.